ROLLY
890.01

Regolamento (UE)
2017/745

IT

WWW.KONG.IT

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INDICE

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Download the digital version

TABLE OF CONTENTS

	1 - SYMBOLS AND SUPPORT	4
	1.1 Symbols	4
	1.2 Support	4
	2 - GENERAL INFORMATION	5
	3 - TECHNICAL FEATURES	6
	3.1 Terminology and materials of the parts	6
	3.2 Technical and performance specifications	8
	3.3 Optional components and spare parts	9
	4 - SPECIFIC INFORMATION	10
	4.1 Intended use	10
	4.2 Preparation	10
	4.3 Setup	11
	4.4 Blocking the patient	16
	4.5 Transporting the patient	18
	5 - MAINTENANCE AND REPAIRS	20
	5.1 General	20
	5.2 Maintenance	20
	5.3 Repair	20
	6 - STORAGE	21
	7 - CHECKS, INSPECTIONS AND SERVICING	22
	7.1 Checks	22
	7.2 Inspections	22
	7.3 Servicing	22
	8 - DURATION	23
	8.1 Duration	23
	8.2 Disposal	23
	8.3 Law obligations	23
	9 - LABELLING AND SYMBOLS	24
	9.1 Medical device label	24
	9.2 Symbols	24
	10 - DOCUMENTS	25
	10.1 Inspections and servicing register	25
	10.2 Maintenance and repair register	26
	10.3 EU Declaration of conformity	27
	11 - LIST OF APPLIED STANDARDS AND REFERENCE	28
	11.1 Applied standards	28
	11.2 Standards used as reference	29

CHAPTER

1

SYMBOLS AND SUPPORT

1.1 SYMBOLS

To make reading this manual comfortable and clear, the following are the symbols used to manage important warnings, for correct and safe use of the device.

REQUISITE FOR CORRECT USE

It identifies the presence of information required to use the device correctly.

INFORMATIVE REQUISITE

It identifies the presence of useful and general information, to guide the user to consciously use the device and/or carry out the actions.

It identifies that the Medical Device is manufactured, designed and produced in compliance with the provisions of the General Safety and Performance Requirements of Regulation (EU) 2017/745 (Class I medical device, in compliance with classification rule 1 as indicated in the attachment VIII).

1.2 SUPPORT

For information, contact Kong Customer Service by:

- telephone +39 0341 630506,
- fax +39 0341 641550,
- email: safetycare@kong.it,

or write to KONG S.p.A. – Via XXV Aprile, 4 – 23804 Monte Marenzo LC - ITALY.

To facilitate assistance operations, always communicate or indicate the serial number (SN) indicated on the label applied to the Medical Device.

CHAPTER 2

GENERAL INFORMATION

The information provided by the manufacturer (hereinafter referred to as information) must be read and clearly understood by the user prior to using the Medical Device (MD hereafter). The information regards the description of the features, performance, assembly, disassembly, maintenance, preservation, disinfection, etc. of the device. Even though the information offers tips on use, this information shall not be considered as a user manual under actual conditions of use.

2.1 WARNINGS AND LIMITATIONS OF USE:

- This device shall be strictly used by people who are physically fit, trained (instructed and taught) to use the device and with specific experience regarding moving the patient or, during training activities, by people under direct supervision of the trainers/ supervisors who guarantee the safety thereof,
- do not use the device before fully reading and understanding this user manual,
- strictly follow the manufacturer's information, improper use of the device is hazardous,
- modifying and/or repairing the device is strictly forbidden,
- all checks described in chapter 7 must be carried out prior to and after using the device. In case of any doubt on the efficiency of the device, the user must replace it immediately,
- non-compliant use, deformations, falls, wear, chemical contamination, exposing textile/ plastic components/devices to temperature below -30°C or above +50°C and metal components/devices to temperatures exceeding 100°C are some examples of causes that can reduce, limit and end the life of the device,
- Before use, ensure that the weight does not exceed the load capacity tested for this MD,
- in order to reduce risks of exposure to / transmission of infectious diseases, clean and disinfect the device as indicated in chapter 5,
- improper use of the patient immobilisation systems can jeopardise the safety of the patient,
- always check the compatibility of the devices used alongside the device by consulting the relative manufacturer's information,
- use of spare parts or optional components different from the ones defined for this MD,
- do not expose the device to sources of heat and at contact with chemical substances. Reduce exposure to direct sunlight as much as possible. At low temperatures and in the presence of humidity, formation of ice could reduce flexibility and increase cutting and abrasion-related risks on textile and synthetic devices,
- report to the manufacturer and the competent authority of the member state where the user and / or patient is established any serious incident occurring in relation to the MD.

All our devices are tested/ checked piece by piece in compliance with the procedures laid down by the Quality System certified in accordance with the UNI EN ISO 9001 standard. Laboratory tests, testing, information and standards do not always reproduce the practical result. Thus, the results obtained under the actual conditions of use of the device may differ, even considerably at times. The best information lies in the continuous practical use under the supervision of skilled/expert/qualified people.

Warning: not suitable for use in an ATEX environment (Directive 2014/34/UE)

CHAPTER 3

TECHNICAL FEATURES

3.1 TERMINOLOGY AND MATERIALS OF THE PARTS

- A - Support surface in high-density polyethylene
- B1 - Slings in polyester
- B2 - Slings in polyester
- C - Slings in polyester with eyelets
- D - Slings in polyester
- E - Sling in polyester
- F - Handles in polyester
- G - Hanging kit in high-tenacity polyester (G1 short part - G2 long part)
- H - Static rope
- I - Buckle with hook in aluminium alloy
- L - Buckle with hook in aluminium alloy
- M - Slings with eyelets
- N - Sling with ratchet
- O - Transport bag.

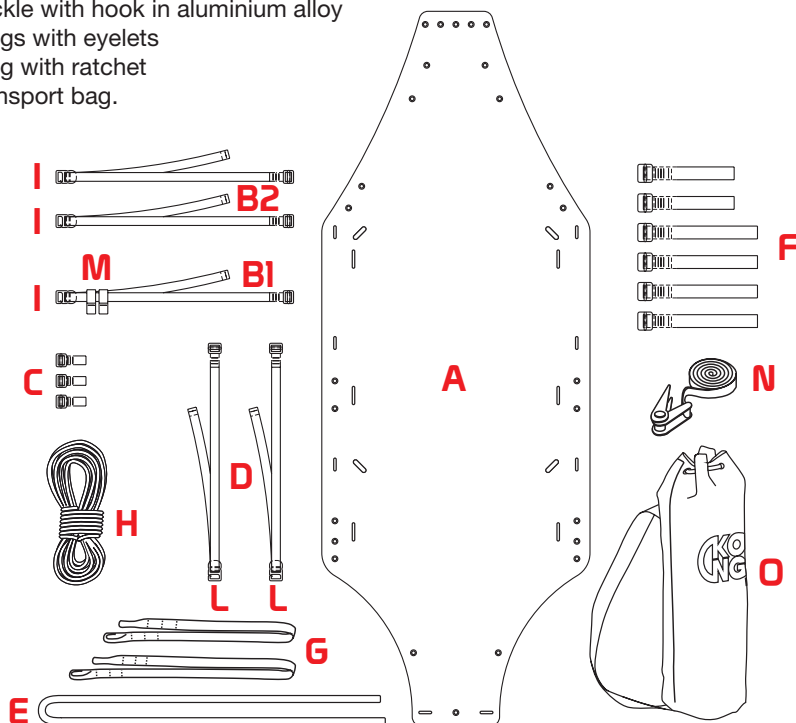
Fig.1

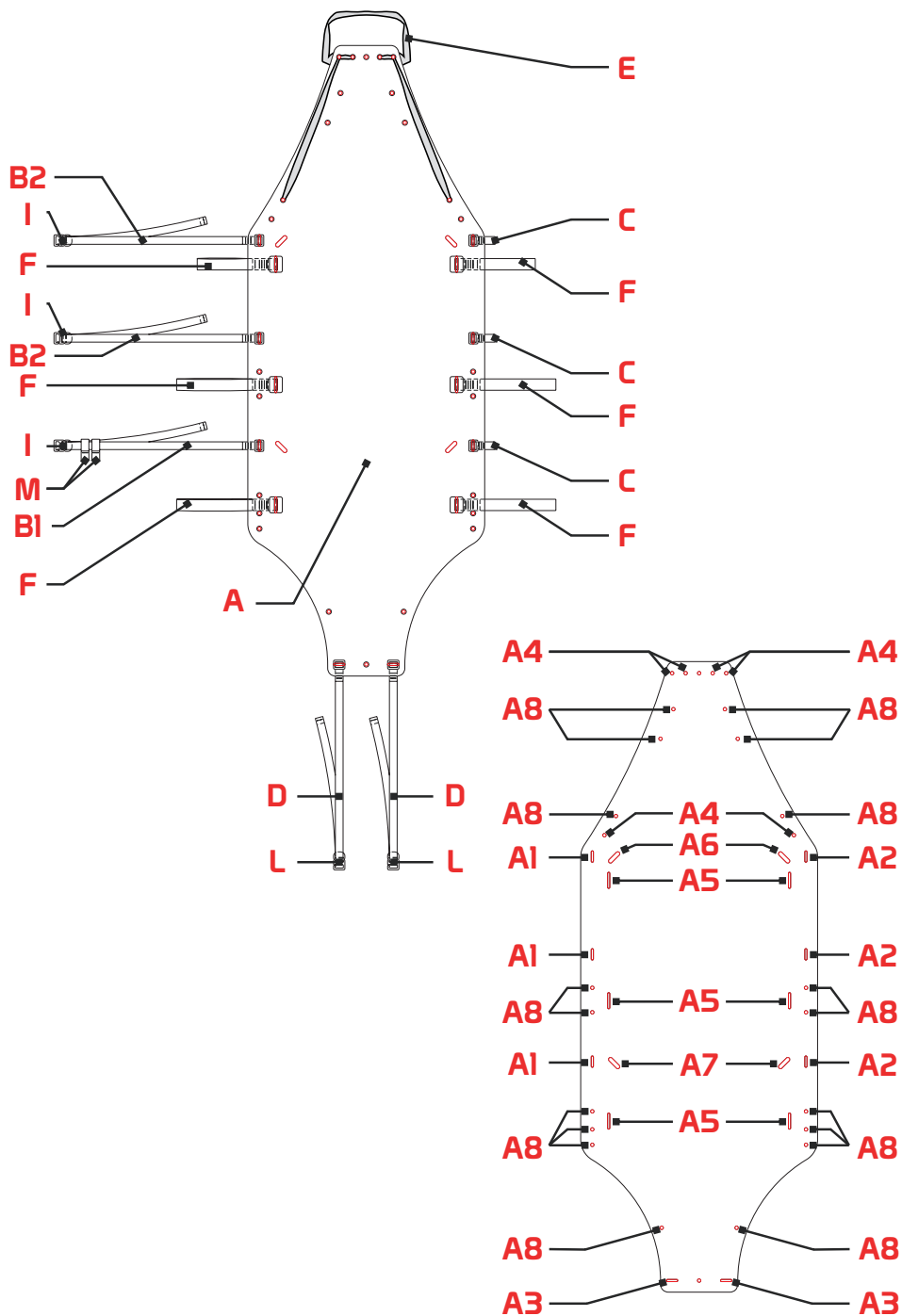
Fig.2

3.2 TECHNICAL AND PERFORMANCE SPECIFICATIONS

3.2.1 Technical Data

Dimensions:	245x92 cm
Dimensions of the DM when stored in the carrying bag:	Ø30 x 110 cm
Total weight of the parts (fig. 1):	9 kg
Rated load:	150 kg

3.2.2 Strength Tests

The MD underwent the following tests to verify the resistance and the conformity to the standard. The results confirm the structural integrity of the MD at the following conditions:

Static resistance Verified by lifting using the side handles (fig. 16):	250 kg
Static resistance Verified by lifting with the Hanging Kit (fig.17):	1500 kg
Static resistance Verified by lifting with the static rope (fig. 18):	1200 kg
Buoyancy force Verified by using the optional kit 880020300KK FLOK-IT (fig. 19):	400 N (~40 kg)

3.3 OPTIONAL COMPONENTS AND SPARE PARTS

The methods of connection and use of the optional components are defined in the relative user instructions.

3.3.1 Optional components

880020300KK FLOK IT

Floating kit for stretchers, made of three cylinders of different lengths (2 side pieces of 1 m each and 1 piece head side of 0,4 m), each of one equipped with two fastening straps with quick release plastick buckles.

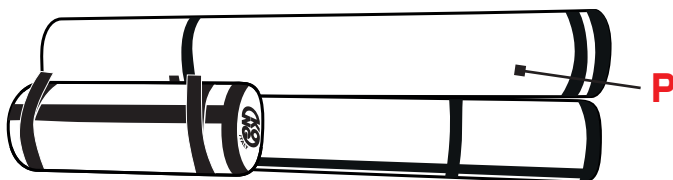
Fig.3

3.3.2 Spare Parts

- B1 - Slings in polyester
- B2 - Slings in polyester
- C - Slings in polyester with eyelets
- D - Slings in polyester
- E - Sling in polyester
- F - Handles in polyester
- G - Hanging kit in polyester
high tenacity
(G1 short part - G2 long part)
- H - Static rope
- N - Sling with ratchet
- O - Transport bag.

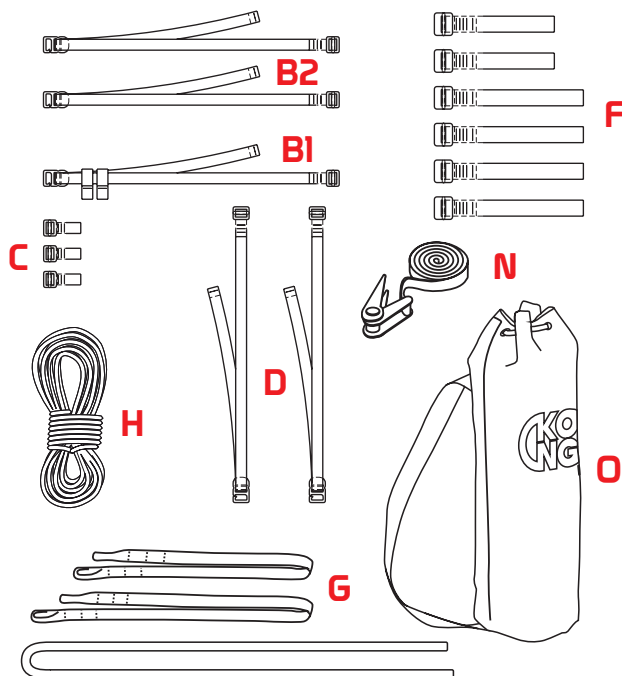
Fig.4

CHAPTER 4

SPECIFIC INFORMATION

4.1 INTENDED USE

The Medical Device (MD) **"ROLLY"** is a rescue sheet that is particularly suited to recovering and transporting a patient - even immobilised on **"X-TRIM"** spinal boards or on the **"VACUUM"** mattresses, in a quick manner and in particularly difficult situations.

Decisions concerning moving and immobilising the patient, as well as the duration, method to be used and combination with other devices, shall be taken and executed by expert and trained personnel only.

The supplied hanging kit (G) makes it possible to use the **"ROLLY"** sheet with a winch.

4.2 PREPARATION

- Remove the "ROLLY" rescue sheet from the transport bag (O) and remove the ratchet sling (N)
- Distend the support surface (A) and flatten it, rolling the ends over in the opposite direction and folding the central section (fig. 5)

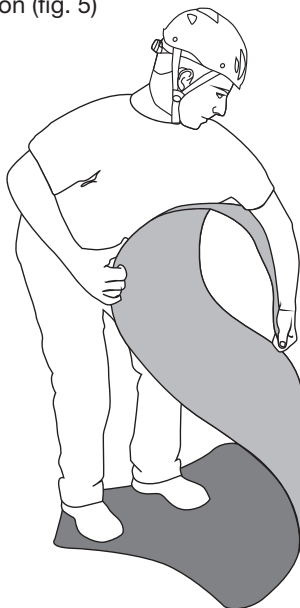
Fig. 5

4.3 SETUP

4.3.1 Setting up for blocking the patient

Prepare the following setup to block the patient:

- Insert the sling (E) in the A4 holes in the support surface (fig. 6).
- Insert the slings (B1 and B2) in the A1 eyelets in the support surface (fig. 7).
- Insert the slings with eyelets (C) in the A2 eyelets in the support surface (fig. 8).
- Insert the slings (D) in the A3 eyelets in the support surface (fig. 9).

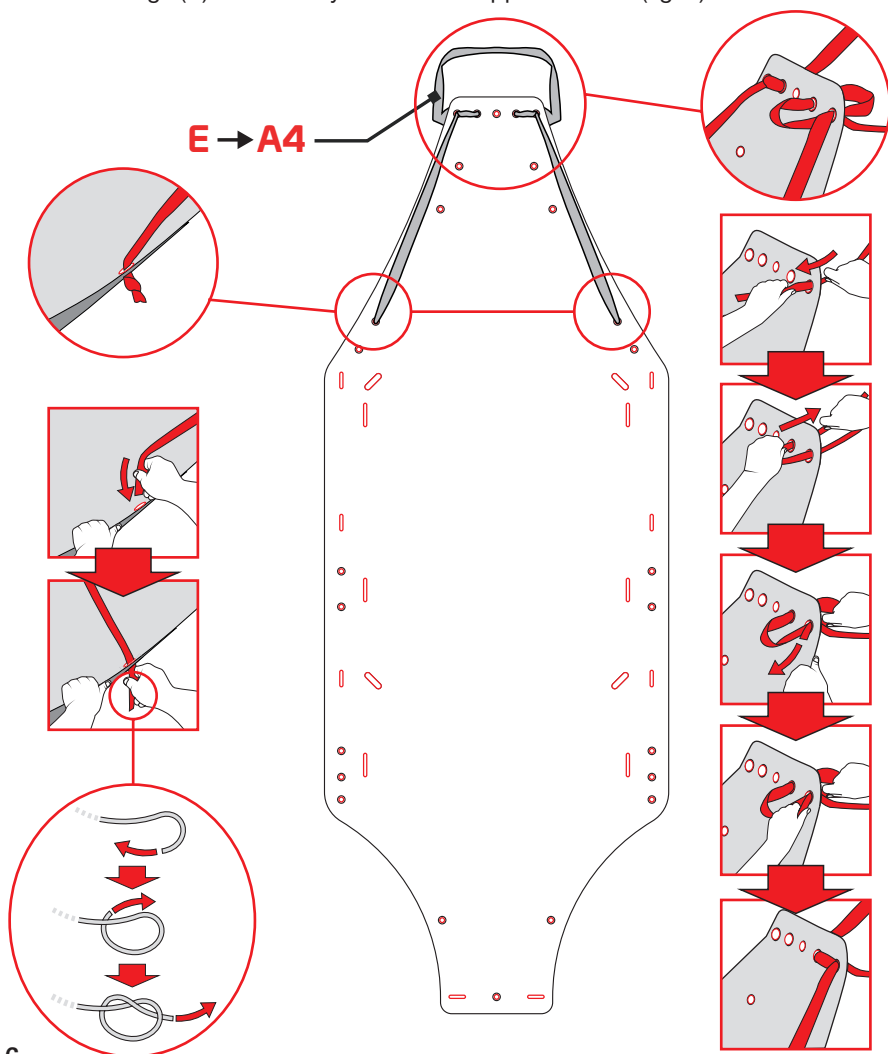
Fig.6

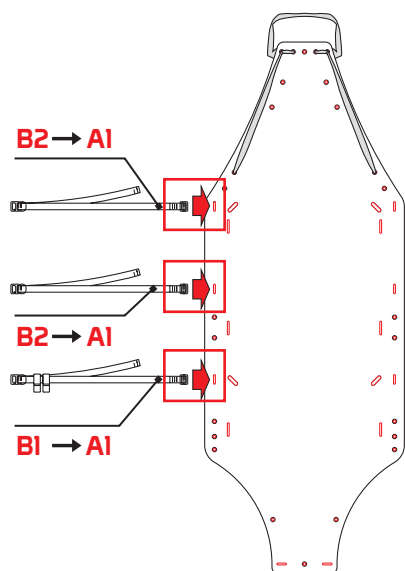
Fig.7

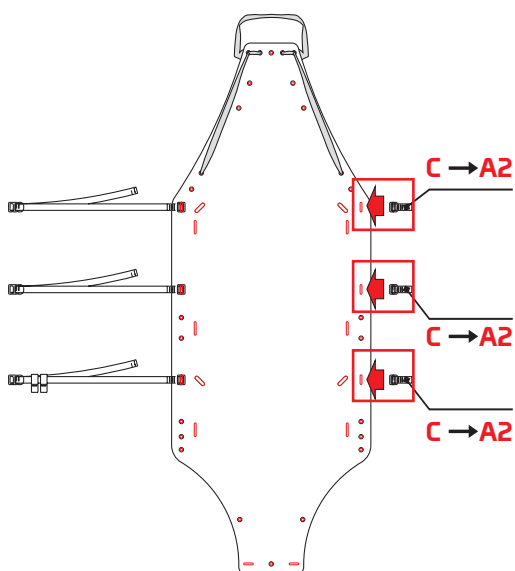
Fig.8


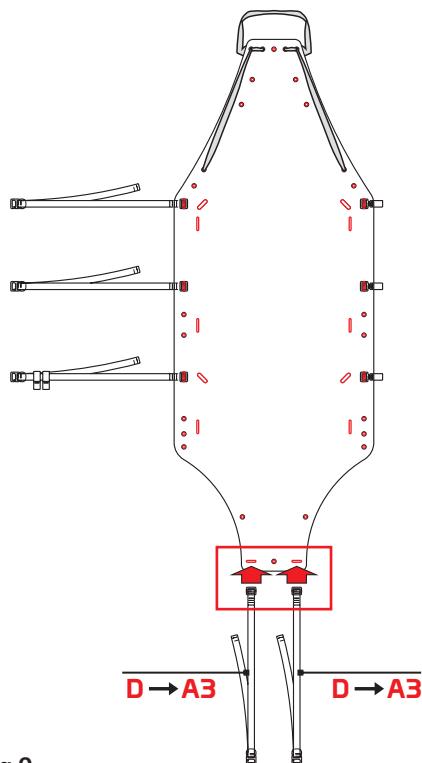
Fig.9

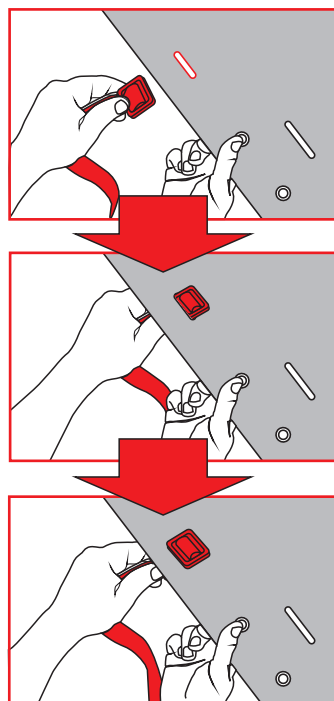4.3.2 Setting up for moving the patient

In relation to the lifting and transport procedures, which the rescue workers considers suitable, set the stretcher up with at least one of the following parts:

a) manual transport: insert only the side handles (F) in the eyelets A5 on the support surface (fig. 10) - Please note: we suggest inserting the shortest couple of handles in the place near the torso of the patient, in order to keep the head higher than the feet.

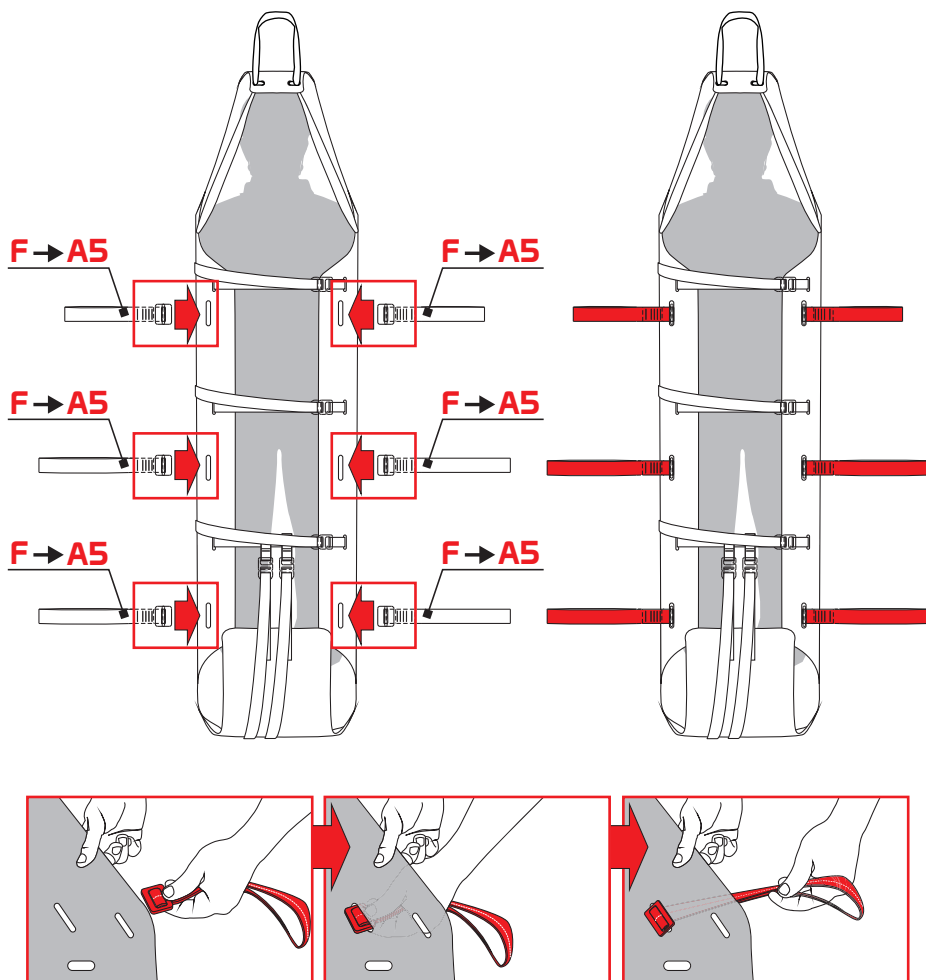
Fig.10

b) lifting/lowering and/or transportation using the lifting kit (G): insert the short hanging kit (G1) in the eyelets A6 and the long hanging kit (G2) in the eyelets A7 on the support surface (fig. 11).

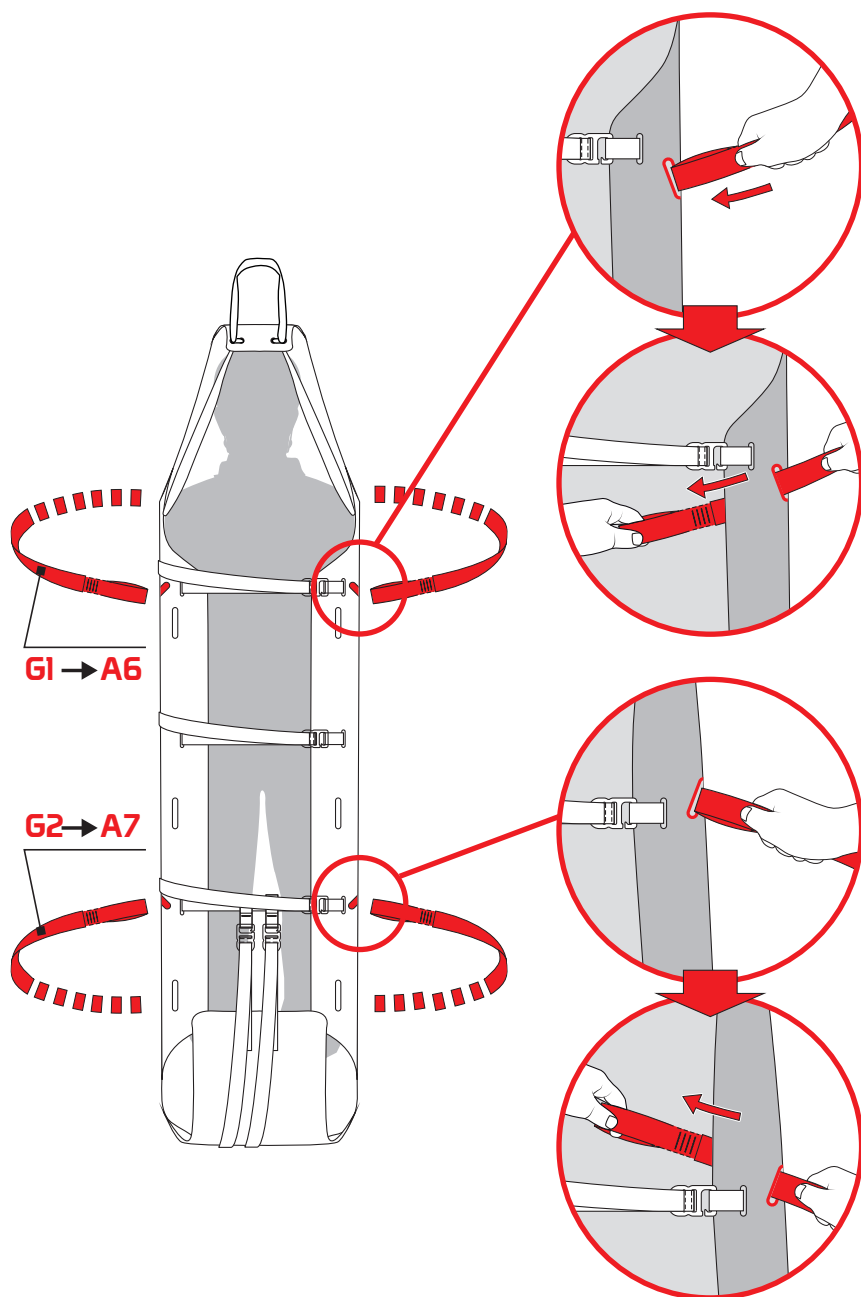
Fig.11

c) lifting/lowering or vertical lowering: insert the rope (H) in the holes A8 in the support surface (fig. 12).

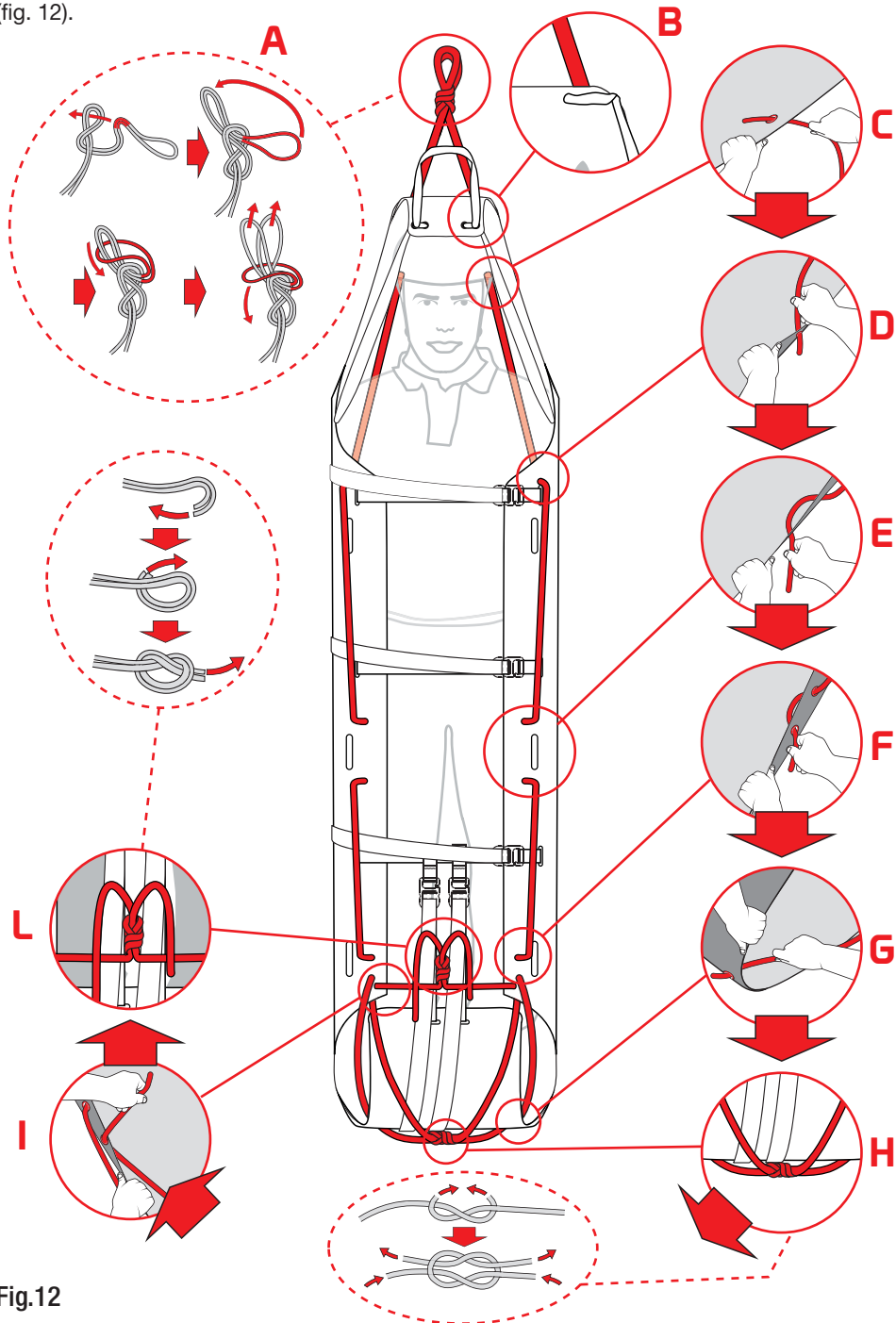
Fig.12

4.4 BLOCKING THE PATIENT

After having positioned the patient on the support surface (A):

- Insert the buckles with hook (I) in the corresponding slings with eyelet (C) and tension the slings (B) to lift the side parts of the support surface (A) enough to contain the patient (fig. 13).
- Insert the buckles with hook (L) in the corresponding slings with eyelet (M) and tension the slings (D) to lift the end part of the support surface (A) enough to contain the patient's feet (fig. 14).


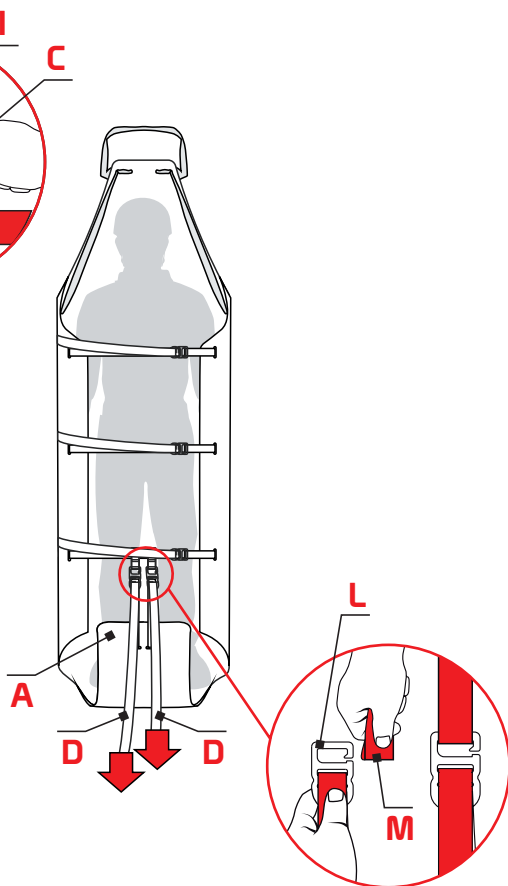
Fig.13

Fig.14

- c) Tension the sling (E) to lift the end part of the support surface (A) enough to contain the head of the patient and immobilise it (fig. 15).

Avoid prolonged contact of the patient's skin with the fabric of the device.

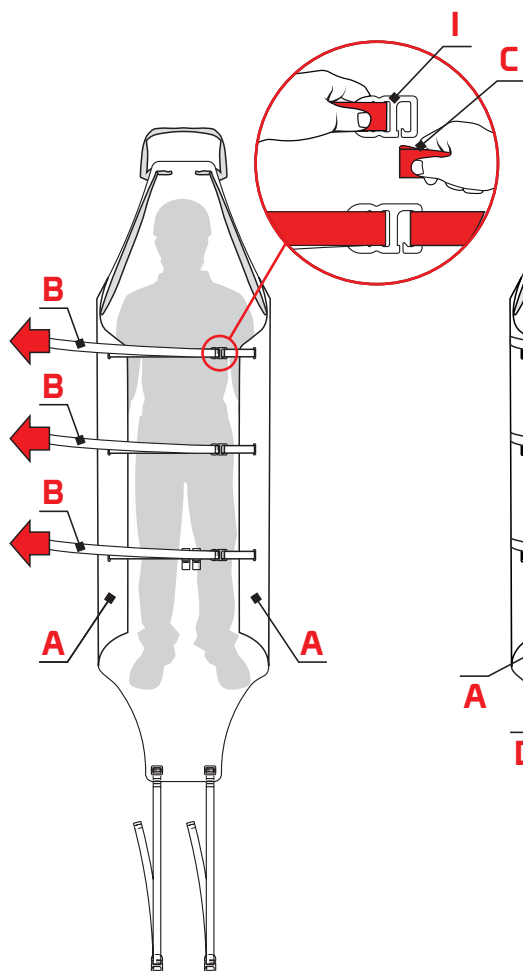This procedure must also be carried out with patients immobilised on the spinal board "X-TRIM" or on the vacuum-operated mattress "VACUUM".

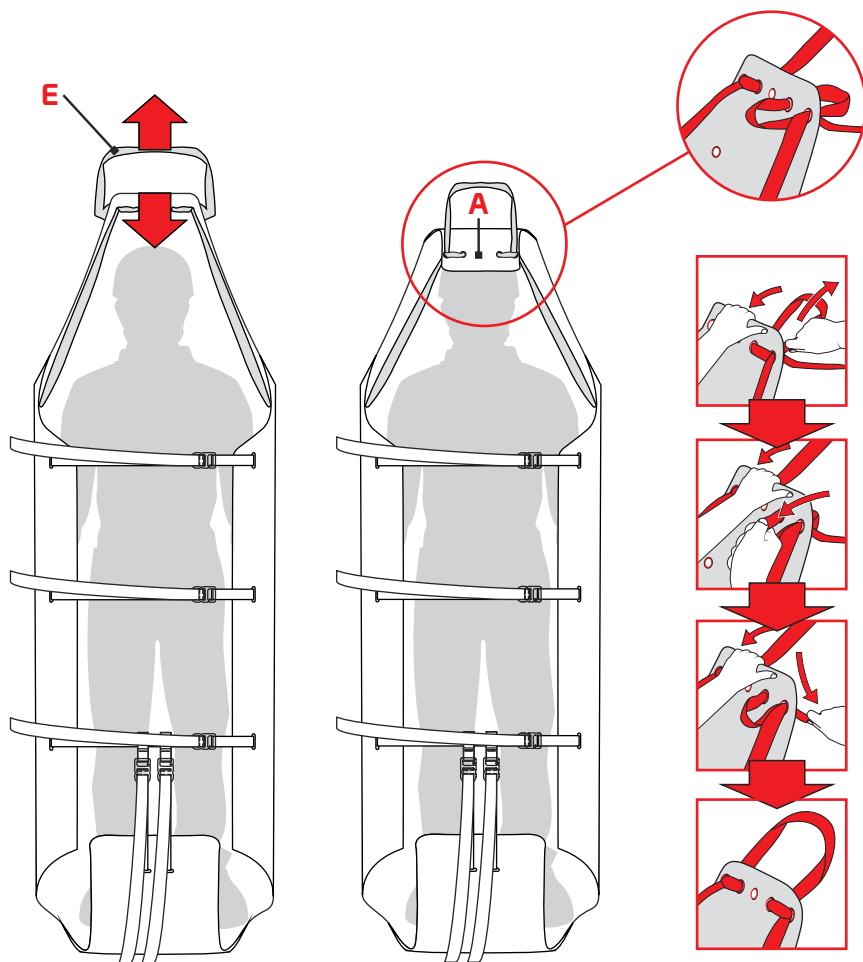
Fig.15

4.5 TRANSPORTING THE PATIENT

The MD “**ROLLY**” is suitable for use in the following modes of use:

- **Manual transportation (fig. 16)** with the handles (F)
- **Lifting and lowering (fig. 17)** with the hanging kit (G)
- **Vertical descent and recovery (fig. 18)** with the rope (H)
- **Water transport (fig. 19)** with the kit 880020300KK FLOK-IT (P)

Before any recovery operation, make sure that the weight does not exceed the tested load capacity!

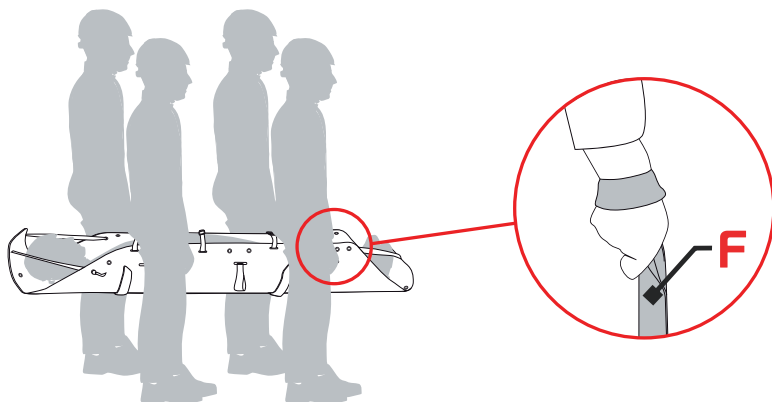
Fig.16

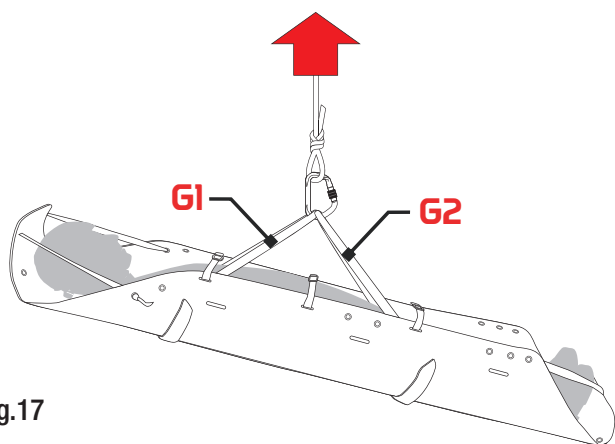
Fig.17

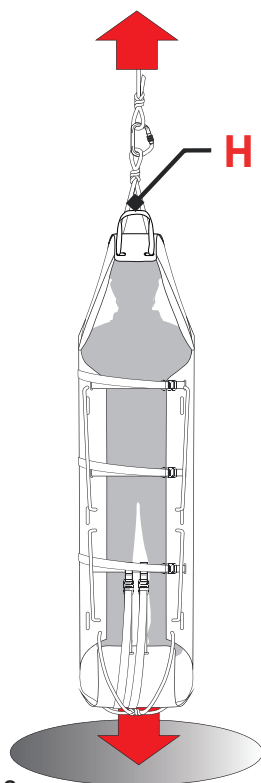
Fig.18

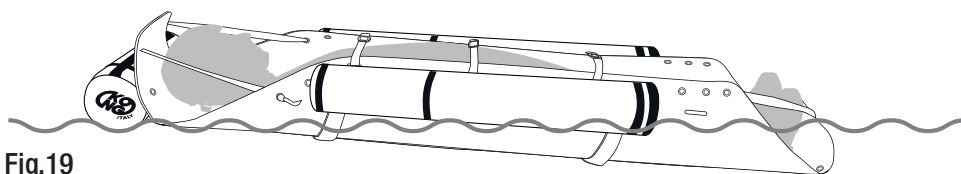
Fig.19

CHAPTER 5

MAINTENANCE AND REPAIR

5.1 GENERAL INFORMATION

The MD is realized with high resistance material against wear and external agents. Despite this, the conditions of use make necessary to perform maintenance interventions, and in certain case, even reparation interventions.

Maintenance and reparation must be recorded on a designated register, an example of which is shown in chapter 10.

5.2 MAINTENANCE

Cleaning: wash the MD by hand in lukewarm water (30-32° C) with soap (Ph 5,5-8,5), then let it dry in the open air and avoid the direct exposure to the sun or source of heat.

do not bleach, do not wash in the washing machine, do not dry in the drying machine, do not stir, do not iron, do not sterilize devices / the textile parts in an autoclave.

Disinfection: submerge the MD for at least 30 minutes in hot water (58-60°C), then let it dry in the open air avoiding the direct exposure to the sun or source of heat.

Do not utilize this method for the disinfection of the product in High Density Polyethylene (E.G. Dyneema®)

5.3 REPARATION

Reparations must be carried out by the manufacturer only. The user is only authorized to change the parts mentioned in the paragraph 3.3.2 with new and original parts.

CHAPTER 6

6 STORAGE

After cleaning, disinfecting (if necessary) and drying, store the MD and its accessories in a dry (relative humidity 40-60%), cool (temperature 10-25°C), dark, well-ventilated place, away from sharp edges, chemicals and other harmful conditions.

Do not store this device wet!

CHAPTER 7

CHECKS, INSPECTIONS AND REVISIONS

7.1 CHECKS

Before and after every use, check the device and ensure that:

- it is appropriate for the use it was intended for,
- it has not been subjected to mechanical deformations and it neither has cracks nor shows signs of wear,
- the textile parts do not have cuts, burns, chemical products residues, excessive hair, wear, in particular check areas at contact with the metal components (buckles, rings, etc.);
- the markings, including labels, are legible.

7.2 INSPECTIONS

The device must be inspected at least once a year, starting from the date of first use, by expert personnel trained and approved by KONG S.p.A. The date of the first use and the outcome of the inspection must be recorded in the Inspections and Servicing Register, regarding which an example is indicated in chapter 10.

7.3 REVISIONS

The device must be serviced by KONG S.p.A., or by expert personnel trained and approved by KONG S.p.A. if:

- any malfunctions are detected,
- the result of the pre and post use check or of the inspections is negative.

The outcome of the servicing must be recorded in the Inspections and Servicing Register, regarding which an example is indicated in chapter 10.

The revised device is guaranteed for one year from the revision date.

CHAPTER 8

DURATION

8.1 DURATION

Carefully read the “Warnings and limitations of use” paragraph in chapter 2.

The duration of this device depends on the outcome of the at least one-yearly inspection (paragraph 7.2) and possible servicing.

The hanging kit, even if not used, must be replaced after 10 years from the first use and in any case no later than 12 years from the date of manufacture.

Important: any devices that have not passed the pre-use, post-use and periodic inspection must be eliminated and made unusable.

8.2 DISPOSAL

Follow the rules in force in the country of use and the waste disposal procedures of the reference hospital facility for proper disposal.


8.3 LAW OBLIGATIONS

Professional and free time activities are often regulated by special national laws which can impose limits and/or obligations concerning the use of these devices. The user has the obligation to know and apply such laws which could provide for limits different from those indicated in this information.

CHAPTER 9

LABELLING AND SYMBOLS

9.1 MEDICAL DEVICE LABEL

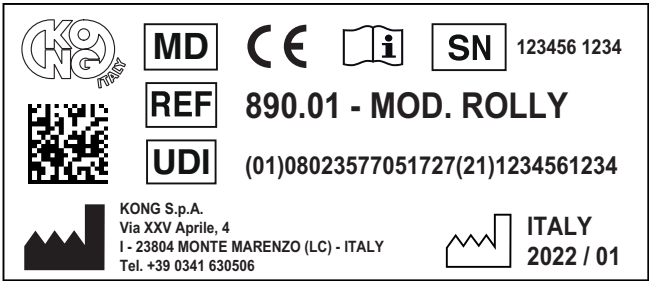Any changes to the positioning of the symbols does not change the content thereof.

9.2 SYMBOLS

Identification and traceability

	Medical Device
	Product identification code
	Unique Device Identifier
	Manufacturer's identification
	Unique serial number
	Country, year and month of production.

Data matrix UDI
Readable with the Apps:
Qbar - Code Scan - UpCode

Information and safety

	Consult the user manual
	Compliant with Regulation (EU) 2017/745

CHAPTER 10

DOCUMENTS

10.1 INSPECTIONS AND SERVICING REGISTER

RESCUE SHEET - Class I MD				
REF.	890.01		MOD.	ROLLY
UDI-DI	08023577051727		SN	
Date of first use				
INSPECTION AND SERVICING				
DATE	I/S	Description	Outcome	Supervisor

10.2 MAINTENANCE AND REPAIR REGISTER

RESCUE SHEET - Class I MD				
REF.	890.01		MOD.	ROLLY
UDI-DI	08023577051727		SN	
Date of first use				
MAINTENANCE AND REPAIRS				
DATE	M/R	Description	Outcome	Supervisor

10.3 EU DECLARATION OF CONFORMITY

The EU Declaration of Conformity is available for download at www.kong.it

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www.kong.it

CHAPTER 11

LIST OF APPLIED AND REFERENCE STANDARDS

11.1 APPLIED STANDARDS

- EN 597-1: 2015
Furniture - Assessment of the ignitability of mattresses and upholstered bed bases - Part 1: Ignition source smouldering cigarette
- EN 1865-1: 2010 + A1: 2015
Patient handling equipment used in road ambulances - Part 1: General stretcher systems and patient handling equipment
- EN 62366-1: 2015 + A1: 2020
Medical devices - Application of usability engineering to medical devices
- EN ISO 846:2019
Plastics - Evaluation of the action of microorganisms
- EN ISO 10993-1: 2025
Biological evaluation of medical devices - Part 1: Evaluation and testing within a risk management process
- EN ISO 10993-5: 2009
Biological evaluation of medical devices - Part 5: Tests for in vitro cytotoxicity
- EN ISO 10993-23: 2021
Biological evaluation of medical devices - Part 23: Tests for irritation
- EN ISO 14971: 2022
Medical devices - Application of risk management to medical devices
- EN ISO 15223-1: 2025
Medical devices — Symbols to be used with information to be supplied by the manufacturer — Part 1: General requirements
- EN ISO 20417: 2021
Medical devices - Information to be supplied by the manufacturer

11.2 STANDARDS USED AS REFERENCE

- EN ISO 13485: 2016
Medical devices - Quality management systems - Requirements for regulatory purposes
- EN 13718-1: 2014 + A1: 2020
Medical vehicles and their equipment - Air ambulances - Part 1: Requirements for medical devices used in air ambulances.

KONG s.p.a.

Via XXV Aprile, 4 - (zona industriale)

I - 23804 MONTE MARENZO (LC) - ITALY

Tel +39 0341630506 - Fax +39 0341641550

www.kong.it